

# Casa Elite

Nursing and residential care bed solutions simplified

Designed with resident safety and comfort in mind

Compliant with  
EN 60601-1,  
EN 60601-1-11\* and  
EN 60601-2-52

Casa Elite  
Standard

Casa Elite  
Low

**CASAELITE**

Designed and developed for the nursing and residential care environment; the Casa Elite Care Home beds provide safety features for both residents and carers combined with an attractive appearance. This range of four section electrical profiling beds with adjustable head, foot and knee sections support effective resident care. The Casa Elite offers multiple profiling and positioning benefits which can reduce other potential complications of immobility.

\*Type B applied parts compliant

**drive**  
**DeVilbiss**  
HEALTHCARE



The Casa Elite is available as a standard or a low model; both bed variants allow the bed to be raised/lowered to a height such that a carer's posture is closer to neutral for the moving and handling of patients.




**CASAELITE**Low

The Casa Elite Low enables the bed platform to be lowered closer to the floor reducing the risk of a fall from height for vulnerable residents. It can also support independent transfers (ingress/egress from the bed) for some individuals.

A large range of side rails also provide additional safety for the resident whilst in bed; with wooden and metal options available including extended rails for the taller resident and height extended rails for those assessed as requiring them.

Available in three wood finishes to compliment room décor. Wood colours available are Beech (this is the standard colour), Light Oak and Walnut.

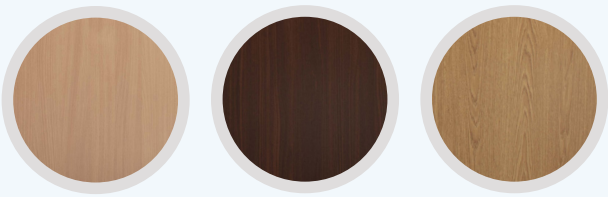

**Beech**  
Standard Colour  
Code: **BECH**

**Walnut**  
Optional Colour  
Code: **WALN**

**Light Oak**  
Optional Colour  
Code: **LOAK**



The Casa Elite comes as standard with an intuitive and easy-to-use 9 button handset.



An optional 10 button handset for Trendelenburg positioning (head down tilt) is available, allowing a tilt of 11 degrees if prescribed for an individual's clinical needs.



The lightweight handset can be locked-out to ensure vulnerable residents' safety.

A simple modular ordering process; choose a standard or low bed, choose your wood colour then add side rails and accessories required.

Interchangeable spare parts and accessories including side rails, grab handles and padding systems, shared across other Drive DeVilbiss Bed Ranges increases flexibility, offers cost and stock holding benefits; whilst reducing training requirements due to familiarity with the shared items, simplifying bed fleet management and maintenance.

**Infusion (Drip) Pole**



**Lifting Pole**



A vast range of compatible accessories are available to meet key user requirements; including a bed extension for the taller residents and both wooden and metal side rails.

The patented full length dipped metal side rails, designed with a stepped down section for easy resident egress and ingress; preventing any excess pressure on the back of the legs whilst moving in or out of the bed.

**Full Length Wooden Side Rails**



**Patented Dipped Side Rails**



Casa Elite Community bed variant



The Casa Elite is supplied on a transport stand to assist with the transportation and storage of the bed.

Casa Elite Community bed variant

The bed is designed for single person installation - all bed sections are below 25kg for safe manual handling. Simple, quick and easy to install requiring only an Allen key (supplied and attached to the bed frame).

The Casa Elite range of beds are compliant to:

**EN 60601-1-11\***  
Medical electrical equipment - Part 1-11 - Requirements for medical electrical equipment and medical electrical systems used in the home healthcare environment.

**EN 60601-1**  
Medical electrical equipment - Part 1 - General requirements for basic safety and essential performance.

**EN 60601-2-52**  
Medical electrical equipment - Part 2-52 - Particular requirements for basic safety and essential performance of medical beds.

\*Type B applied parts compliant



Designed by our bed experts in Halifax, West Yorkshire; the Casa Elite bed range provides many clinical, safety and comfort benefits to ensure optimal support to residents and carers.

Flexibility and choice

Wood finishes:

- 1. The standard wood finish for the Casa Elite is Beech - alternative colours are Light Oak or Walnut. This allows the wood finish to match the room décor.

Side rails:

- 2. The bed can be specified at the point of ordering or retro-fitted later to accommodate changing clinical needs with a range of side rails. This includes 2 bar wooden side rails or 2 bar dipped or standard metal side rails.

The patented metal dipped side rails are innovative as they are designed for easy patient egress and ingress; minimising excess pressure on the backs of the legs whilst transferring into or out of the bed.

- 3. Interchangeable accessories including side rails, grab handles and padding systems, shared across other Drive DeVilbiss bed ranges. Increasing flexibility, whilst reducing training requirements due to familiarity with the shared items.

- 4. The Casa Elite can be specified with a ratchet to adjust the leg/calf section angle. The calf section of the bed can be raised to enable the users lower legs to remain elevated. Once the knee break has been set at the desired angle, the leg stay angle can be set by lifting and releasing this section manually on the ratchet. To release, raise the leg section to the highest height and then lower.



Auto-Contour



Reverse Trendelenburg



Trendelenburg



Cardiac chair position



- 5. Reverse Trendelenburg function (foot down position) places the body on an incline with the head elevated up to 11°, profiling the patient to maintain a socially interactive and comfortable position in the bed.

- 6. The Cardiac chair position additionally provides a comfortable seated position for users nursed in bed, further promoting interaction and comfort whilst undertaking activities, by utilising the auto contour and foot down tilt in combination.

- 7. Auto contour function is a pre-programmed profiling position that raises/lowers both the backrest and knee break simultaneously to promote comfort and minimise migration in the bed. *SOURCE; Preece J 'Total bed management: The way forward in pressure sore prevention', British Journal of Nursing 1999 8 (22) 1524-1529*

- 8. Intuitive and easy to use 9 button handset is provided as standard with the Casa Elite. The handset includes a lock-out facility to prevent unauthorised use for additional patient safety.



An optional 10 button handset for Trendelenburg positioning (head down tilt) is available, allowing a tilt of 11° if prescribed for an individuals clinical needs.

- 9. Low height option available at just 19.2cm minimum height from the floor to the mattress platform; reducing the risk of a fall from height for vulnerable residents

- 10. Optional bed platform extension accommodates taller individuals by increasing the length by 16cm. Extended side rails are also available (wooden or metal) to be used in conjunction with the platform extension to support and provide safe, comfortable and dignified care.



# How to order your Casa Elite Bed

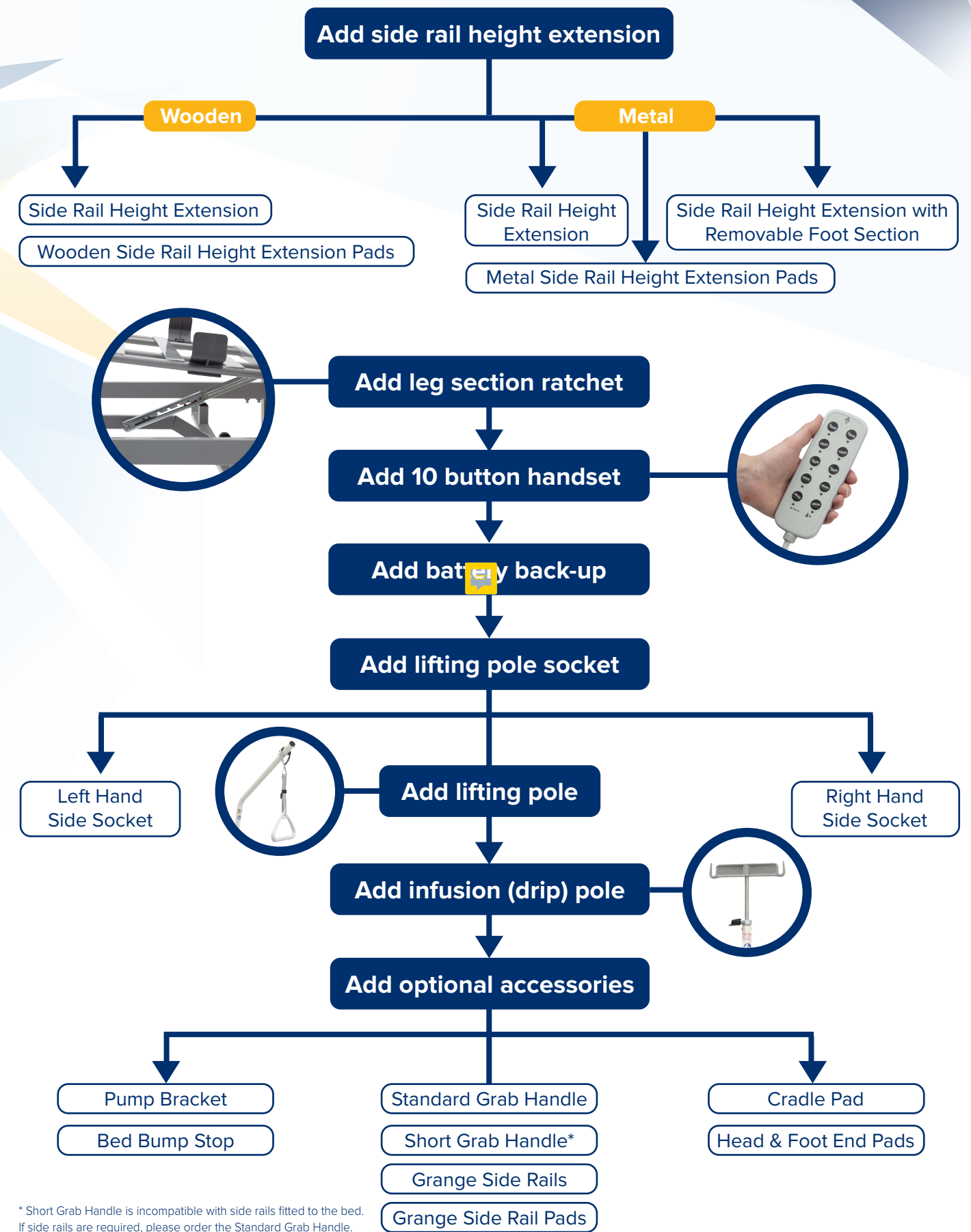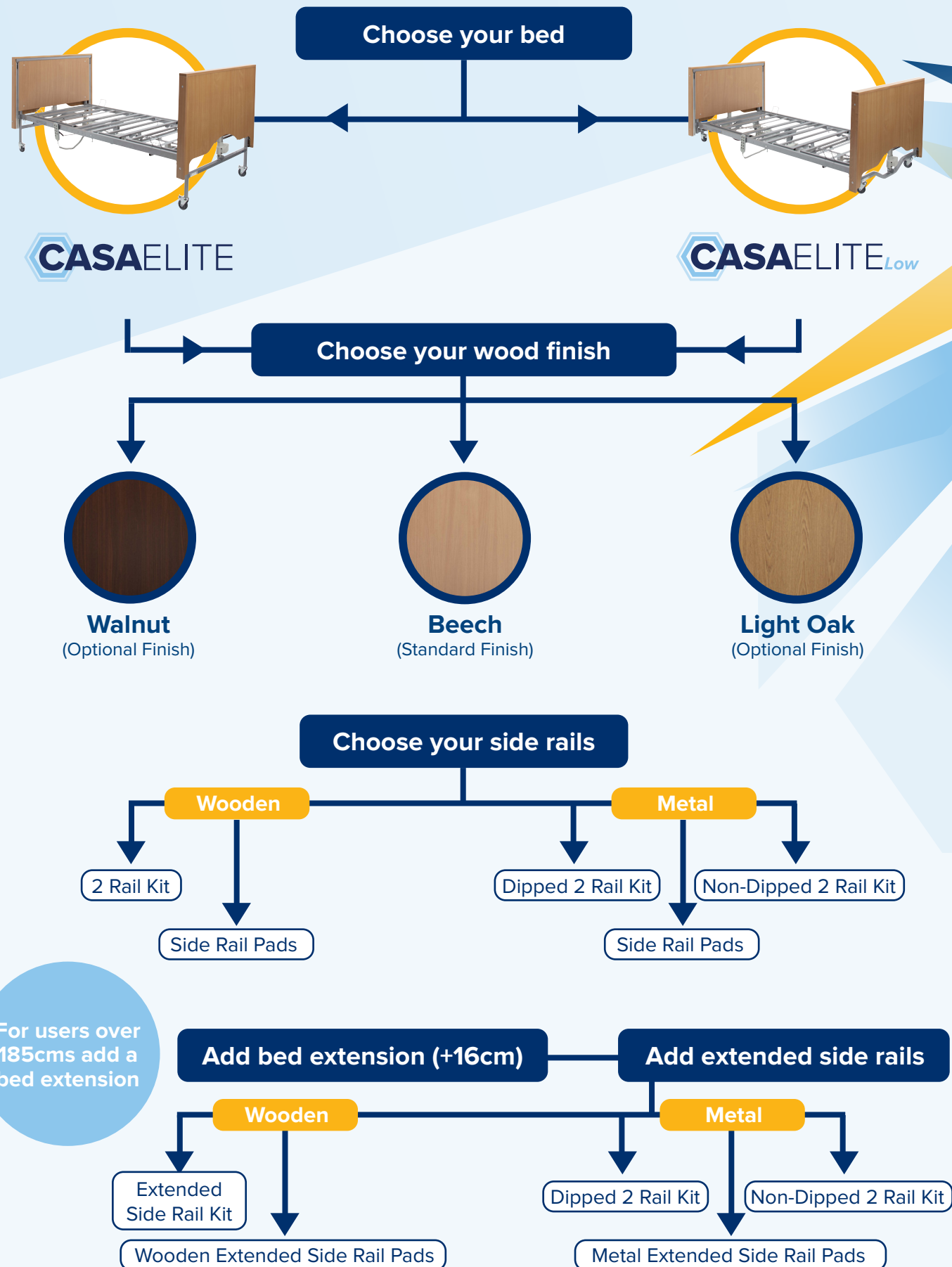

## Specifications

\*Add the wood colour to the end of the bed product code (Beech = **BECH**, Light Oak = **LOAK** and Walnut = **WALN**)

	
CASA/C/STD/*COLOUR*	CASA/C/LOW/*COLOUR*
Overall Bed Width: <b>102.5cm</b>	Overall Bed Width: <b>102.5cm</b>
Overall Bed Length: <b>214cm</b>	Overall Bed Length: <b>214cm</b>
Overall Platform Width: <b>90.7cm</b>	Overall Platform Width: <b>90.7cm</b>
Overall Platform Length: <b>200.2cm</b>	Overall Platform Length: <b>200.2cm</b>
Mattress Platform Height (min): <b>38.6cm</b>	Mattress Platform Height (min): <b>19.2cm</b>
Mattress Platform Height (max): <b>80.6cm</b>	Mattress Platform Height (max): <b>61.2cm</b>
Underbed Clearance for Hoist: <b>34.2cm</b>	Underbed Clearance for Hoist: <b>15cm</b>
Trendelenburg Tilt <b>11° +/- 1°</b>	Trendelenburg Tilt <b>11° +/- 1°</b>
Max User Weight: <b>165kg (26st)</b>	Max User Weight: <b>165kg (26st)</b>
Max Safe Working Load: <b>200kg (31½st)</b>	Max Safe Working Load: <b>200kg (31½st)</b>
Product Weight on Transport Stand: <b>66kg</b>	Product Weight on Transport Stand: <b>66kg</b>
Warranty: <b>1 year full parts and labour + 2 years parts only</b>	Warranty: <b>1 year full parts and labour + 2 years parts only</b>

## Kitted Codes

\*Add the wood colour to the end of the bed product code (Beech = **BECH**, Light Oak = **LOAK** and Walnut = **WALN**)

Product Code	Description
CASA/K/C/STD/*COLOUR*	Casa Elite Care Home Beds Standard with Wooden Side Rail Kit
CASA/K/C/LOW/*COLOUR*	Casa Elite Care Home Beds Low with Wooden Side Rail Kit

## Metal Side Rails

Product Code	Description
RAIL/M/KIT/D/S	Metal Dipped Side Rail Kit including Finger Kit - Silver
RAIL/M/KIT/D/EX/S	Metal Extended Dipped Side Rail Kit including Finger Kit - Silver
RAIL/M/KIT/ND/S	Metal Non-Dipped Side Rail Kit including Finger Kit - Silver
RAIL/M/KIT/ND/EX/S	Metal Extended Non-Dipped Side Rail Kit including Finger Kit - Silver
RAIL/M/HEX/S	Metal Side Rail Height Extension - Silver
RAIL/M/HEX/RF/S	Metal Side Rail Height Extension with Removable Foot Section - Silver

Leaflet Code: **LL489**

Drive DeVilbiss Healthcare Ltd.  
Sidhill Business Park, Holmfield, Halifax, West Yorkshire HX2 9TN

Casa\_Elite\_Care\_Home\_Bed\_Brochure\_LL489\_220210\_2.0

## Wooden Side Rails

\*For wooden side rails add the chosen wood colour to the end of the bed product code (Beech = **BECH**, Light Oak = **LOAK** and Walnut = **WALN**)

Product Code	Description
RAIL/W/KIT/*COLOUR*	Wooden Side Rail Kit including Finger Kits
RAIL/W/KIT/EX/*COLOUR*	Wooden Extended Side Rail Kit including Finger Kits
RAIL/W/HE/EX	Wooden Side Rail Height Extension

## Pads

Product Code	Description
1275/PS	Wooden Side Rail Pads
1275/PS/EXT	Wooden Extended Side Rail Pads
SOLITE/PRO/SR/PAD	Metal Side Rail Pads
SOLITE/PRO/SR/EXT/PAD	Metal Extended Side Rail Pads
SOLITE/PRO/SR/HE/1/PAD	Integral Metal Side Rail Height Extension Pads
1316/PS/CREAM	Grange Side Rail Pad
SOLITE/PRO/EPS/PAD	Cradle Pad
BED/HF/PAD	Bed Foot and Head Board Pads

## Accessories

Product Code	Description
CASA/BED/EX	Casa Bed Length Extension 16cm
BED/RATCHET	Leg Section Ratchet
<del>CASA/BATTERY</del>	<del>Casa Battery Back-up Facility</del>
CASA/PUMPBKT	Casa Pump Bracket
BED/LP	Lifting Pole
CASA/LPSOCKET/RH	Lifting Pole Socket, Right Hand
CASA/LPSOCKET/LH	Lifting Pole Socket, Left Hand
BED/DP	Infusion (Drip) Pole
BED/BUMPSTOP	Bed Bump Stop
1210/GRAB	Standard Grab Handle
1216/GRAB**	Short Grab Rail
GRANGE/SR	Grange Side Rail

\*\* 1216/GRAB is incompatible with side rails fitted to the bed. If side rails are required, please order 1210/GRAB.